

Centra Line™

Central Venous Catheter Sets

Biometrix offers a broad range of catheter kits and trays for central venous catheterization. We provide safe and reliable devices for a wide variety of clinical interventions. Centra-Line™ enables several procedures to take place simultaneously via a single puncture site: CVP monitoring, continuous and simultaneous infusion, blood sampling, blood product administration and TPN administration.



Customized CVC Packs

- customized drapes
- surgical gowns
- gauze
- swabs
- hypodermic needles
- needle counters

Adult Sets

4 Lumen Catheter Sets | Adult

Code	Description	Needle	Neutron®	Quick Insertion Syringe
HH-8515	8.5Fr x 15cm	Super Sharp	-	-
HH-8520	8.5Fr x 20cm	Super Sharp	-	-
HH-8530	8.5Fr x 30cm	Super Sharp	-	-
HH-I8515	8.5Fr x 15cm	Y-Needle	-	-
HH-I8520	8.5Fr x 20cm	Y-Needle	-	-
HH-I8530	8.5Fr x 30cm	Y-Needle	-	-
HH-E8515	8.5Fr x 15cm	Y-Needle	+	-
HH-E8520	8.5Fr x 20cm	Y-Needle	+	-
HH-E8530	8.5Fr x 30cm	Y-Needle	+	-
HH-R8516	8.5Fr x 15cm	Super Sharp Echogenic	+	+
HH-R8521	8.5Fr x 20cm	Super Sharp Echogenic	+	+
HH-R8531	8.5Fr x 30cm	Super Sharp Echogenic	+	+

All sets contain a Nitinol Kink Resistant Guidewire

3 Lumen Catheter Sets | Adult

Code	Description	Needle	Neutron®	Quick Insertion Syringe
HJ-7015	7Fr x 15cm	Super Sharp	-	-
HJ-7020	7Fr x 20cm	Super Sharp	-	-
HJ-7030	7Fr x 30cm	Super Sharp	-	-
HJ-I7015	7Fr x 15cm	Y-Needle	-	-
HJ-I7020	7Fr x 20cm	Y-Needle	-	-
HJ-I7030	7Fr x 30cm	Y-Needle	-	-
HJ-E7015	7Fr x 15cm	Y-Needle	+	-
HJ-E7020	7Fr x 20cm	Y-Needle	+	-
HJ-E7030	7Fr x 30cm	Y-Needle	+	-
HJ-R7016	7Fr x 15cm	Super Sharp Echogenic	+	+
HJ-R7021	7Fr x 20cm	Super Sharp Echogenic	+	+
HJ-R7031	7Fr x 30cm	Super Sharp Echogenic	+	+
HJ-R8516	8.5Fr x 15cm	Super Sharp Echogenic	+	+
HJ-R8521	8.5Fr x 20cm	Super Sharp Echogenic	+	+

All sets contain a Nitinol Kink Resistant Guidewire

2 Lumen Catheter Sets | Adult

Code	Description	Needle	Neutron®	Quick Insertion Syringe
HK-7015	7Fr x 15cm	Super Sharp	-	-
HK-7020	7Fr x 20cm	Super Sharp	-	-
HK-7030	7Fr x 30cm	Super Sharp	-	-
HK-I7015	7Fr x 15cm	Y-Needle	-	-
HK-I7020	7Fr x 20cm	Y-Needle	-	-
HK-I7030	7Fr x 30cm	Y-Needle	-	-
HK-E7015	7Fr x 15cm	Y-Needle	+	-
HK-E7020	7Fr x 20cm	Y-Needle	+	-
HK-E7030	7Fr x 30cm	Y-Needle	+	-
HK-R7016	7Fr x 15cm	Super Sharp Echogenic	+	+
HK-R7021	7Fr x 20cm	Super Sharp Echogenic	+	+
HK-R7031	7Fr x 30cm	Super Sharp Echogenic	+	+

All sets contain a Nitinol Kink Resistant Guidewire

1 Lumen Catheter Sets | adult

Code	Description	Needle	Neutron®	Quick Insertion Syringe
H0-4015	14Ga x 15cm	Super Sharp	-	-
H0-4020	14Ga x 20cm	Super Sharp	-	-
H0-4030	14Ga x 30cm	Super Sharp	-	-
H0-6015	16Ga x 15cm	Super Sharp	-	-
H0-6020	16Ga x 20cm	Super Sharp	-	-
H0-6030	16Ga x 30cm	Super Sharp	-	-
H0-8015	18Ga x 15cm	Super Sharp	-	-
H0-8020	18Ga x 20cm	Super Sharp	-	-
H0-I4015	14Ga x 15cm	Y-Needle	-	-
H0-I4020	14Ga x 20cm	Y-Needle	-	-
H0-I4030	14Ga x 30cm	Y-Needle	-	-
H0-I6015	16Ga x 15cm	Y-Needle	-	-
H0-I6020	16Ga x 20cm	Y-Needle	-	-
H0-I6030	16Ga x 30cm	Y-Needle	-	-
H0-E4015	14Ga x 15cm	Y-Needle	+	-
H0-E4020	14Ga x 20cm	Y-Needle	+	-
H0-E4030	14Ga x 30cm	Y-Needle	+	-
H0-E6015	16Ga x 15cm	Y-Needle	+	-
H0-E6020	16Ga x 20cm	Y-Needle	+	-
H0-E6030	16Ga x 30cm	Y-Needle	+	-
H0-R4016	14Ga x 15cm	Super Sharp Echogenic	-	+
H0-R4021	14Ga x 20cm	Super Sharp Echogenic	-	+
H0-R4031	14Ga x 30cm	Super Sharp Echogenic	-	+
H0-R6016	16Ga x 15cm	Super Sharp Echogenic	-	+
H0-R6021	16Ga x 20cm	Super Sharp Echogenic	-	+
H0-R6031	16Ga x 30cm	Super Sharp Echogenic	-	+

All sets contain a Nitinol Kink Resistant Guidewire



Pediatrics Sets



3 Lumen Catheter Sets | Pediatrics

Code	Description	Neutron®	Quick Insertion Syringe
HJ-5008	5.5Fr x 8cm	-	-
HJ-5013	5.5Fr x 13cm	-	-
HJ-5030	5.5Fr x 30cm	-	-
HJ-E5008	5.5Fr x 8cm	+	-
HJ-E5013	5.5Fr x 13cm	+	-
HJ-E5030	5.5Fr x 30cm	+	-
HJ-R5008*	5.5Fr x 8cm	+	+
HJ-R5013*	5.5Fr x 13cm	+	+

* Contains Nitinol Kink Resistant Guidewire

2 Lumen Catheter Sets | Pediatrics

Code	Description	Neutron®	Quick Insertion Syringe
HK-4005	4Fr x 5cm	-	-
HK-4008	4Fr x 8cm	-	-
HK-4013	4Fr x 13cm	-	-
HK-4030	4Fr x 30cm	-	-
HK-5006	5Fr x 6cm	-	-
HK-5008	5Fr x 8cm	-	-
HK-5013	5Fr x 13cm	-	-
HK-5020	5Fr x 20cm	-	-
HK-E4005	4Fr x 5cm	+	-
HK-E4008	4Fr x 8cm	+	-
HK-E4013	4Fr x 13cm	+	-
HK-E5006	5Fr x 6cm	+	-
HK-E5008	5Fr x 8cm	+	-
HK-E5013	5Fr x 13cm	+	-
HK-R4005*	4Fr x 5cm	+	+
HK-R4008*	4Fr x 8cm	+	+
HK-R4013*	4Fr x 13cm	+	+
HK-R5006*	5Fr x 6cm	+	+
HK-R5008*	5Fr x 8cm	+	+
HK-R5013*	5Fr x 13cm	+	+

* Contains Nitinol Kink Resistant Guidewire

1 Lumen Catheter Sets | Pediatrics

Code	Description	Neutron®	Quick Insertion Syringe
HO-9010	20Ga x 10cm	-	-
HO-9013	20Ga x 13cm	-	-
HO-9015	20Ga x 15cm	-	-
HO-9020	20Ga x 20cm	-	-
HO-7008	24Ga x 8cm	-	-
HO-7015	24Ga x 15cm	-	-
HO-7020	24Ga x 20cm	-	-

Catheter Tray Contents

- Radiopaque PU catheter
- Guidewire
- Introducer needle / Y-Introducer needle / Super Sharp Echogenic Needle
- 5ml syringe/Quick insertion syringe
- Vessel dilator
- Heparin locks/Needleless valves
- Scalpel



Quantity per adult box: 20/Quantity per pediatric box: 10




Centra Line™ Central Venous Catheter Sets



Kink Resistant Guidewires

The Biometrix Nitinol (kink resistant) guidewires are designed to withstand sharp angles without buckling and coiling, thus facilitating passage of the CVC through the patient's vascular system.

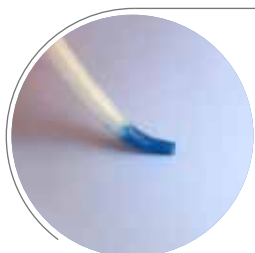

Flex Tip

The geometric tip design along with the high-quality biocompatible material allows a safe and easy insertion and manipulation through the vascular system



Needle Free Valves

Designed to prevent needle stick injuries and reduce infection, the design creates a safe and easy access to the patient. Such a valve can be used with any male luer syringe.



Y-introducer Needle

The Biometrix Y-introducer needle offers a unique design, which allows access to the vein without disconnecting the syringe attached to the needle. The Y-introducer needle has a one-way integral valve for preventing the risk of blood leakage and embolism during insertion.



Quick Insertion Syringe

By using the Biometrix Quick Insertion Syringe the guidewire can be inserted directly into the hollow plunger/barrel, through the introducer needle and into the vascular system. The Quick Insertion Syringe is simpler and faster for the physician, less traumatic for the patient and also serves as a fluid containment device to help protect from infection and cross contamination.